



Small Protein Purification on SRT-10 100Å SEC Prep Size Column

SRT-10 SEC-100 100Å 21.2 x 300mm
(PN: 225100-21230)



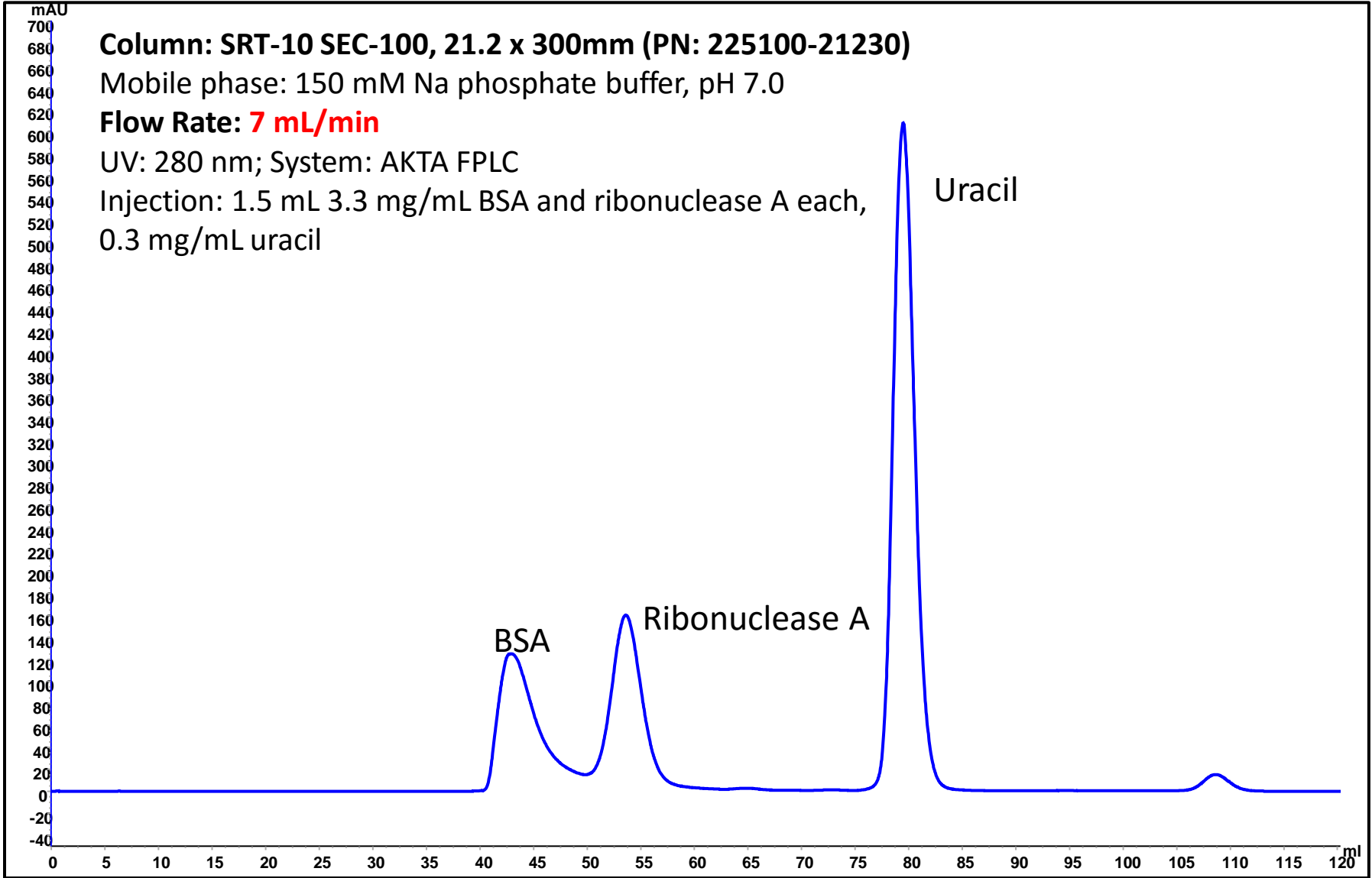

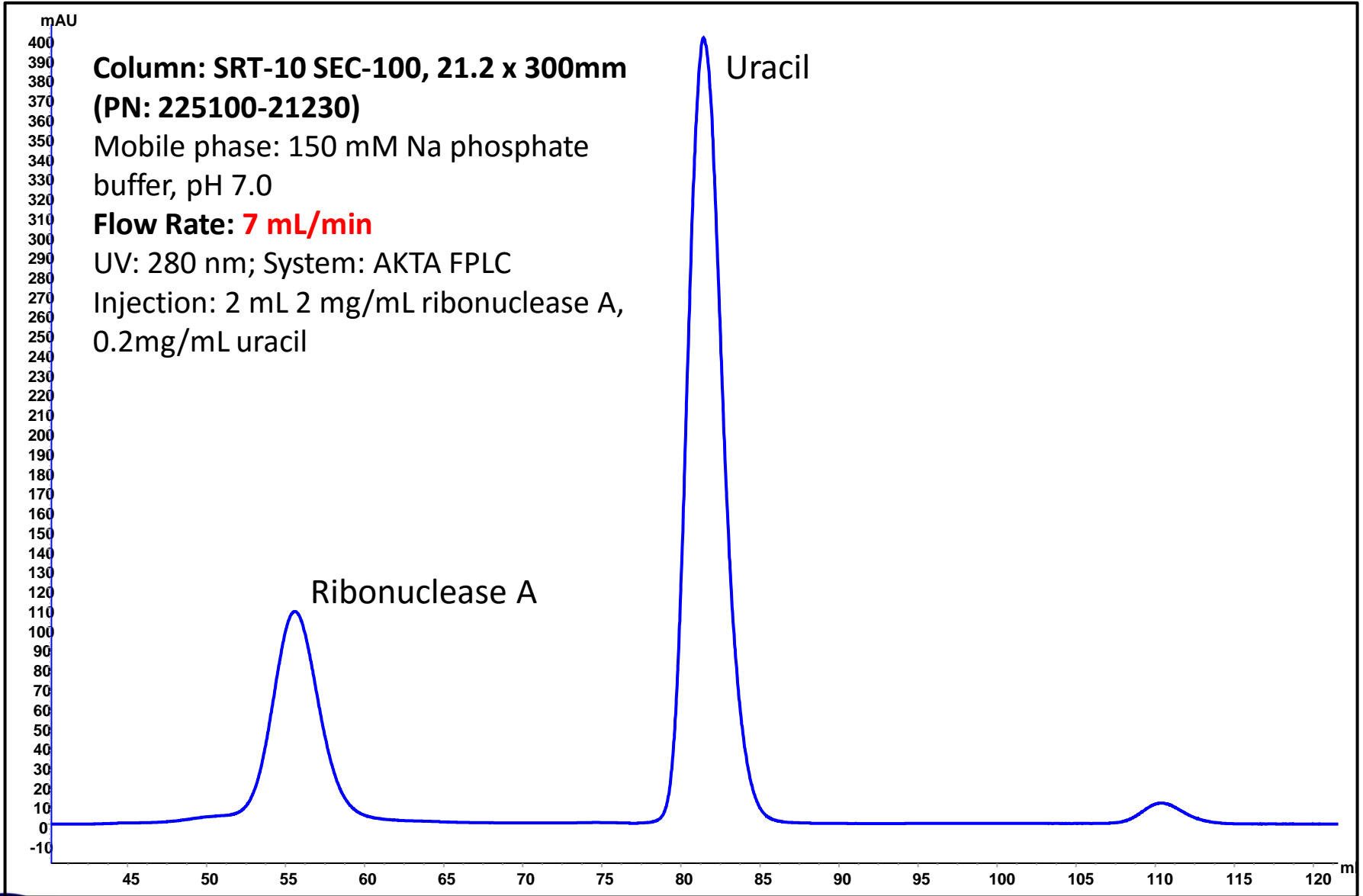

Column: SRT-10 SEC-100, 21.2 x 300mm (PN: 225100-21230)

Mobile phase: 150 mM Na phosphate buffer, pH 7.0

Flow Rate: 7 mL/min

UV: 280 nm; System: AKTA FPLC

Injection: 1.9 mL 2 mg/mL BSA, 0.2 mg/mL uracil

